

HARRIS HEARING AID CENTER

CARE FOR YOUR HEARING AIDS

Clean Your Hearing Aids

- Use a soft, dry cloth or tissue to wipe the entire hearing aid and the ear bud/ear mold. See Fig. 1 and Fig. 2
- Use the cleaning brush to remove any wax that the tissue did not remove. See Fig. 3 Brush the entire hearing aid, including the:

- Ear bud or ear mold
- Push button
- Battery door
- Sound tube attachment
- Microphone
- Receiver








- It is best to **wipe and brush the hearing aid immediately** after you take it out of your ear because the wax is still warm and soft, it has not had time to dry and harden.

Remove and Clean the Sound Tube

- Wax can plug up the sound tube. If the sound tube is plugged the sound will be reduced or **completely blocked off, with no feedback at all**. If the **sound tube is plugged you will get feedback** from the hearing aid **when the sound tube is removed**.
- To remove the sound tube for cleaning:
 - Secure the sides** of the aid with your fingers.
 - Hold the hearing aid firmly while twisting the sound tube.
 - Twist the tube clockwise** a quarter turn. See Fig. 4
 - Pull the tube away** from the hearing aid. See Fig. 5
- After the sound tube is removed: See Fig. 6
 - Insert the tube cleaner** into the end of the sound tube that connects to the hearing aid. **Push the cleaner through** the sound tube **until it can be pulled out the opposite end**.

Store Your Hearing Aids

When not wearing your hearing aids **store them in a cool, dry place**. Open the battery door or **remove the battery**. Keep your hearing aids dry. **Do not expose them to extreme temperatures**. **Children, animals, excessive heat, hairspray, water, moisture and high humidity can damage your hearing aids.**

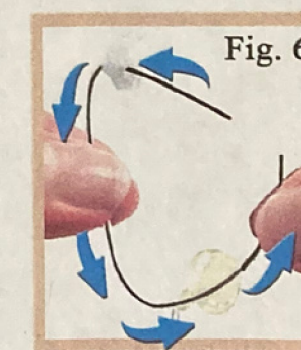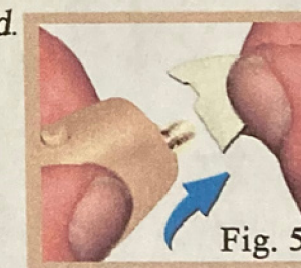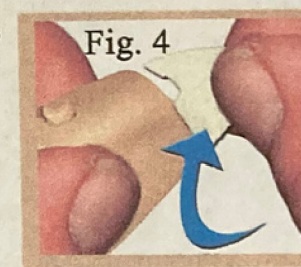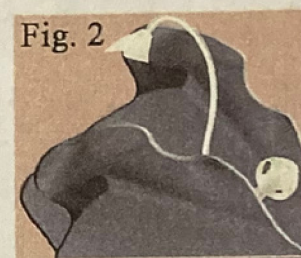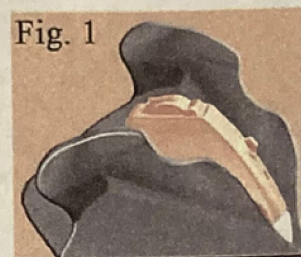

If you have any questions or concerns about hearing aids or instructions please call us at

801-373-6827